


Benzoyl peroxide determination

Summary



Benzoyl peroxide (BPO) is added to bone cement as a radical initiator for the methyl methacrylate monomer. BPO is highly reactive and degrades readily, with a temperature dependent half-life of one hour at 92 °C and one minute at 131 °C. The degradation products, benzoic acid and oxygen, are not significantly toxic, but the BPO remaining to initiate the polymerization is decreased, leaving unreacted monomer in the cement, and the potential for incomplete cure.

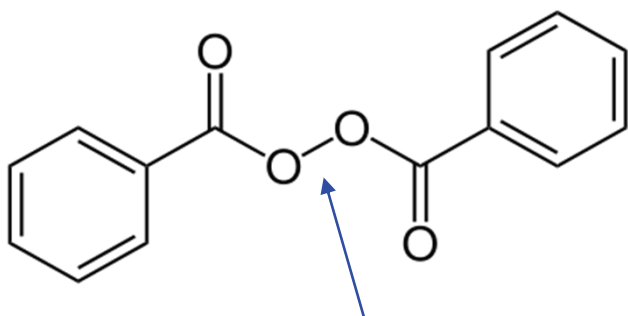

Benzoyl peroxide (BPO) splits at the O-O bond forming two radicals

Description

Common methods for the measurement of BPO in the powder are through HPLC, UV-Vis or titration. A method for titration developed at Cambridge Polymer Group, based on ASTM F451, resulted in data in close agreement with that obtained through HPLC, as shown, but with improved simplicity and substantially reduced cost. Although HPLC will continue to be the gold standard, this titration method allows simple wet-chemistry to be used to quantitatively determine BPO concentration in a sample. Since BPO is known to degrade with time, this assay will prove invaluable for long-term tracking of bone cement efficacy and storage-induced deterioration.

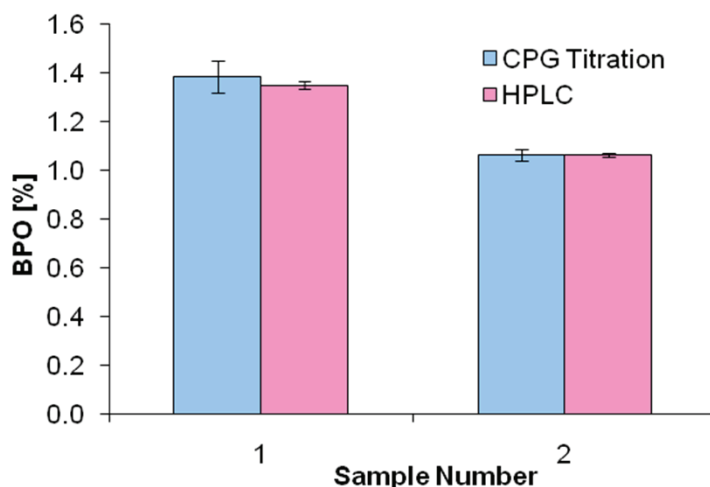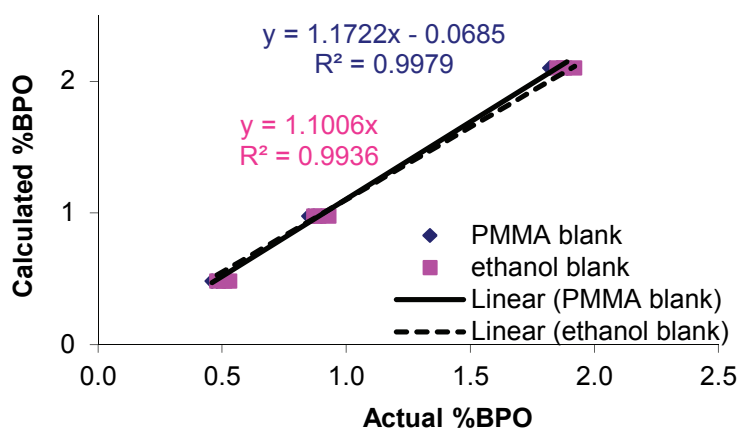

Uses

- Reaction chemistry
- Quality control
- Biomedical materials
- Bone cements



ANALYZE • RESEARCH • CREATE

Cambridge Polymer Group, Inc. is a contract research laboratory specializing in materials. We partner with our clients to solve problems utilizing our multi-disciplinary research team and full service laboratory.

We work with clients throughout the product life cycle to:

- **Develop new materials**
- **Design prototypes for proof-of-concept studies**
- **Create and execute experimental design**
- **Validate and verify manufacturing processes**
- **Perform root-cause analysis in product failures**

Cambridge Polymer Group, Inc. was founded in 1996 to provide a cost-effective resource for testing, research and development to clients who need periodic access to Ph.D.-level scientists and their support structure. We have developed a host of testing methods and materials for our clients, which number more than 300.



Cambridge
Polymer Group



Cambridge
Polymer Group

Cambridge Polymer Group, Inc.
56 Roland Street, Suite 310
Boston, MA 02129

Ph: 1 (617) 629 4400
Fax: 1 (617) 629 9100

info@campoly.com
www.campoly.com
ISO 9001:2008 Certified

Rev: 06/03/11