


Parallel microsample melt rheology

Summary

Conventional methods of determining the flow of polymer melts are expensive and time consuming, and require skilled operators and large sample volumes. In addition, quality control frequently requires multiple specimens for accuracy. CPG has developed a multi-specimen, low-volume, high-throughput melt rheometer for QC and development that will find use with high-value materials or in environments where parallel sampling is desirable.

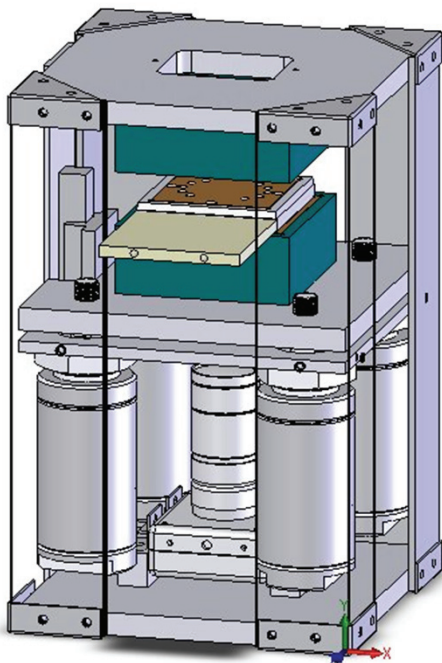


Description

In 2010 NIST published the development of a simpler, pressure drive rheometer that utilize a slit geometry to determine the rheology of a polymer melt. This original concept was extend by CPG with SBIR grant 9.02.04-5.TT to allow enhanced functionality, greater flexibility and improved throughput. A digital camera system is used to monitor the free-surface flow front of a polymer melt as it is driven through a standard slit geometry. This system is theoretically robust and allows multiple samples to be run at high temperatures in a design that can be easily assembled and disassembled. The device should find usage in quality control environments or in high-added-value markets such as combinatorial chemistry and speciality materials.

Specifications

Parameter	Value
Slit height x width	0.132 x 1.3 mm
Slit Length	13 mm
Shear rate	0.001-10,000 1/s
Temperature	-40 – 450 °C
Drive pressure	~0 – 12.4 MPa
Resolution	0.022 mm (0.17% strain)
Sample volume	2.2 μ L
Sample count	5-10

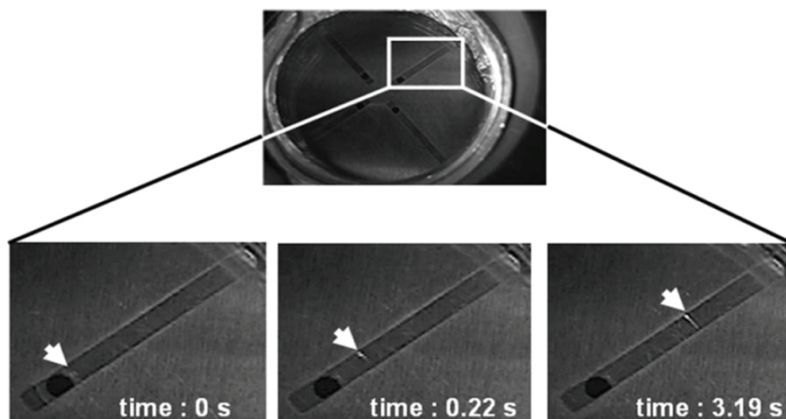

Images: Moon et al. JoR 2008

Markets

- Combinatorial polymer development
- Protein solutions in pharmaceutical development
- Quality control environments



**Cambridge
Polymer Group**

Cambridge Polymer Group, Inc.
56 Roland Street, Suite 310
Boston, MA 02129

Ph: 1 (617) 629 4400
Fax: 1 (617) 629 9100

info@campoly.com
www.campoly.com
ISO 9001:2008 Certified

Rev: 06/03/11

ANALYZE • RESEARCH • CREATE

Cambridge Polymer Group, Inc. is a contract research laboratory specializing in materials. We partner with our clients to solve problems utilizing our multi-disciplinary research team and full service laboratory.

We work with clients throughout the product life cycle to:

- **Develop new materials**
- **Design prototypes for proof-of-concept studies**
- **Create and execute experimental design**
- **Validate and verify manufacturing processes**
- **Perform root-cause analysis in product failures**

Cambridge Polymer Group, Inc. was founded in 1996 to provide a cost-effective resource for testing, research and development to clients who need periodic access to Ph.D.-level scientists and their support structure. We have developed a host of testing methods and materials for our clients, which number more than 300.



Cambridge
Polymer Group



Cambridge
Polymer Group

Cambridge Polymer Group, Inc.
56 Roland Street, Suite 310
Boston, MA 02129

Ph: 1 (617) 629 4400
Fax: 1 (617) 629 9100

info@campoly.com
www.campoly.com
ISO 9001:2008 Certified

Rev: 06/03/11